

aPrehliadky — user guide

3D documentation of bodies, persons, objects and scenes with LiDAR · version 2.2 · forensika.eu

1. What you get and what the app needs from you

aPrehliadky builds a 3D model of a deceased body, a living person, an object or a scene from ordinary iPhone frames. With LiDAR every frame gets depth and gravity, the model is built directly in metres and can be measured. Everything runs on the phone: no cloud, no account, no upload, no AI-completed geometry.

Output: a virtual examination, findings anchored to the surface, measurements (distance, path along the surface, area, volume), surroundings from LiDAR (room mesh), a PDF report, a 3D documentation certificate (class A–D) and an .f3d package with hashes as an attachment for the case file.

Three things decide the quality: 1) the calibration target in the frame, 2) a slow, complete walk around the subject at three heights, 3) even light without flash. Without the target a measurement cannot be independently defended.

The app is a documentation tool. It supports professional judgement and does not replace it. It is not a certified measuring instrument — accuracy depends on capture conditions and is verified with a reference. The user is responsible for the lawfulness of capture and handling of the documentation.

2. Before capture — checklist

- 1 Device: iPhone Pro or iPad Pro with LiDAR, iOS 18+, charged to at least 50 % (processing takes minutes and heats the phone).
- 2 Calibration target T2 printed at 100 % (scan list → menu □ → Printable calibration target T2). Check with callipers: scale 100 mm, QR A 40 mm, fiducial centres B–D 48.0 mm. Without correct printing the reference is worthless. The T2 target has four QR fiducials — the app measures them in 3D from LiDAR and compares six distances with the reference (automatic metrology, written into the report and the certificate).
- 3 Lay the target flat next to the subject, in the plane of the captured surface, not at the edge of the frame. For a body beside the trunk, for an object right next to it.
- 4 Light: even, diffuse (ceiling light, overcast sky). Turn the flash off. Direct sun interferes with LiDAR; hard shadows spoil the texture.
- 5 Surface: dry glossy and wet areas if possible; transparent and uniform surfaces are hard for photogrammetry — LiDAR holds the geometry there.
- 6 Space: walk freely around the subject. If you cannot (body against a wall), the hidden side will be missing in the model — note it in the remarks.
- 7 Sound: the app signals progress by sound — theme Deep Space (default: a quiet deep tone and an airy whistle, unobtrusive) or Star Trek (bridge sensor warble); if it is not appropriate on site, choose No sound with the speaker icon on the capture screen. The pitch rises with coverage; a new sector and completion have their own signal.

3. New scan — what the options mean

- 1 Tap +. Enter the case reference, examiner and place — they go into the report and the package name.
- 2 Profile: Deceased — examination (densest capture, 200 frames, accuracy) · Living person (quick capture, 120 frames, the person must not move) · Object / item (stone, tool, trace; circle at three heights).
- 3 Extent: Overview — whole subject, or Detail — region (wound, tattoo, hand). The triangle budget is the same, so Detail gives a much finer model. Whole body and injury detail are two separate scans.
- 4 Mode: Silent — 4K video (recommended: no shutter sound, depth exactly synchronous with the image) or Photo — full resolution (12 MP frames, slower pace).
- 5 Also save surroundings from LiDAR (on by default): the room mesh is collected during capture — position of the body in the room, distances from walls and objects.
- 6 Press Start capture.

4. Capture step by step

- 1 Aim at the centre of the subject (for a body at the trunk) and tap Set scan centre. Coverage is computed around this point.
- 2 Walk slowly at about 0.8–1.2 m (object 0.3–0.5 m). Frames trigger themselves by translation and rotation of the camera — you press nothing.
- 3 Watch the radar top right: 36 sectors × 3 heights. White wedge = where you stand; coloured = covered. Green radar = coverage is sufficient. The sound rises with coverage, a new sector adds a spark, completion plays a chord.
- 4 Cover three heights: from waist height, from shoulder height and from above. Every point of the surface should appear on at least three frames from different angles.
- 5 If “Moving too fast” or “Few details in view” appears, stop, wait for recovery and continue more slowly.
- 6 Tap Finish capture (at least 20 frames, ideally once the radar is green). Processing starts by itself.

4.1 Deceased body — examination

- Target beside the trunk, visible from all sides. Start at the head, circle the whole body, then a second lap at shoulder height and a third from above (arm stretched over the body).
- Slow down at the head, hands and feet and add close frames (30–40 cm) — the most frequent sites of findings.
- Injuries you know you will measure: capture additionally as a separate scan with extent Detail and the target right next to them.

4.2 Living person

- The person stands or lies still, arms slightly away from the body, gaze fixed. The Living person profile uses fewer frames so capture takes under a minute.
- Chest breathing usually does not matter; limb movement does — repeat the scan if the person moved.

4.3 Object

- Matte, contrasting background (not a glossy table), target next to the object. Circle at three heights: from below, level, from above. Capture small objects close up with extent Detail.
- Document the underside with a second scan after turning the object over.

4.4 Scene / room

- After circling the body, sweep the camera slowly over the walls, floor and nearby objects — the “surroundings ... □” counter on screen grows. This fills the room mesh (chapter 9).
- For processing choose extent Subject with surroundings or Whole scene if the photogrammetric model should include the surroundings.

Most frequent capture mistakes: fast movement (blurred frames), skipped top height (crown, back of the hand), target out of frame or tilted, flash, capturing against a window. Each of them degrades measurement more than any processing setting.


5. Processing

Reconstruction runs entirely on the device, minutes (200 frames 10–20 min). The screen stays on and dimmed; leave the phone on the charger and undisturbed. A notification arrives on completion and the examination opens.

- Reconstruction extent: Subject (crop to the volume around the centre — the whole triangle budget goes into the body) · Subject with surroundings · Whole scene. The whole scene remains available as a second view.
- Sensitivity: High helps with pale skin and low-texture surfaces (takes longer).
- Recompute at any time without new capture — e.g. another extent or sensitivity.
- Processing runs on iPhone and iPad with A17 Pro / M1 or newer. On the Mac processing is not available; a package with a finished model opens on the Mac for examination, measuring and the report.

6. Examination — controls

- One finger rotates, two fingers pan, pinch zooms. View menu: front, back, left, right, top, bottom.
- The bottom panel has two modes: Measure and Mark finding. The switch says what a tap on the surface does. Nothing is hidden in gestures — a long press is only a shortcut to the full finding form.
- The list icon top right opens all findings and measurements; the palette icon turns the texture off (bare geometry shows holes and artefacts); the share icon opens export.
- AR 1 : 1 (ARKit icon, iPhone/iPad only): place the model on a table at real size and compare with a real ruler.

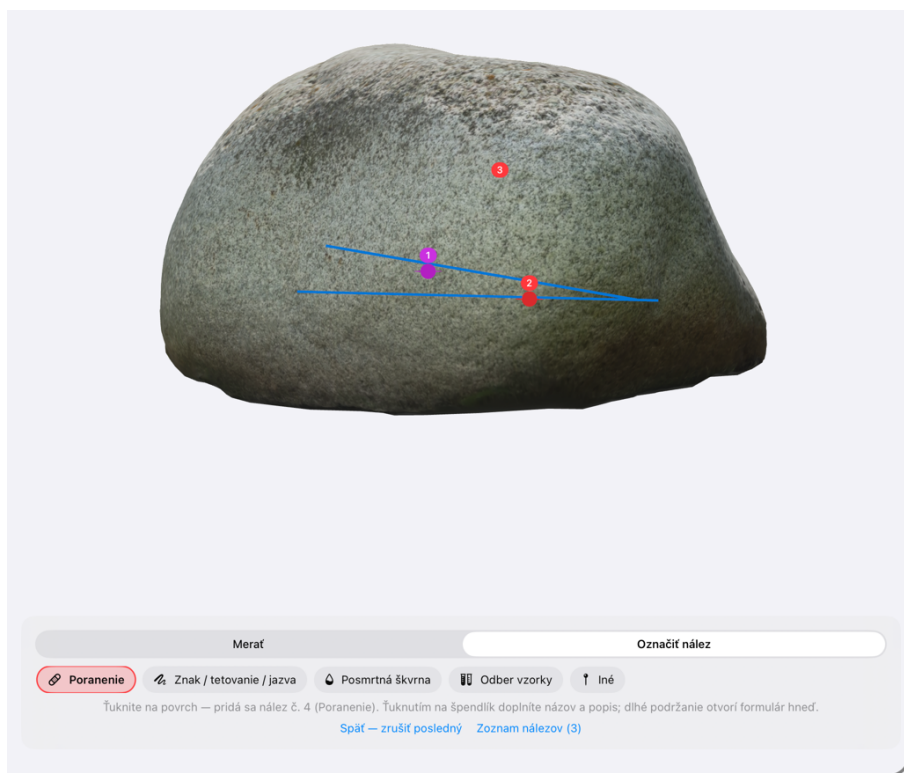


Mark finding mode: choose the type once, every tap on the surface adds a numbered pin (here no. 3).


7. Findings — how to mark precisely and describe

- 1 Switch the panel to Mark finding and choose the type: Injury · Mark / tattoo / scar · Livor mortis · Sample · Other.
- 2 Zoom in so the finding fills half the screen and tap exactly its centre (for a wound the centre of the wound, for a bruise the centre of the area). The pin gets a sequence number and the type.
- 3 Tap the pin → Title (short, as in the autopsy report: “laceration right forehead”) and Note (location relative to landmarks, dimensions from measurement, colour, margins, depth, direction).
- 4 Mistake: Undo — remove the last one deletes the last pin. A misplaced older finding: tap it → Move on the model → tap the correct place.
- 5 Numbering is preserved; in the report findings are listed in order with their position in the model (x, y, z in metres), so they can be found again.

Mark while zoomed in, never from afar. One finding = one pin at its centre; dimensions belong to measurement (chapter 8), not to an estimate. Title = what and where, note = how it looks. Use the same words as in the autopsy report — a colleague must pair them without explanation.

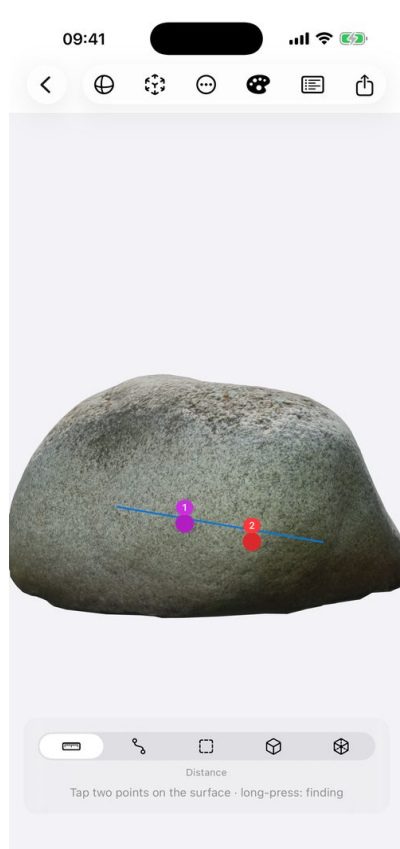

Stone — examination with two findings (labels 1 and 2)

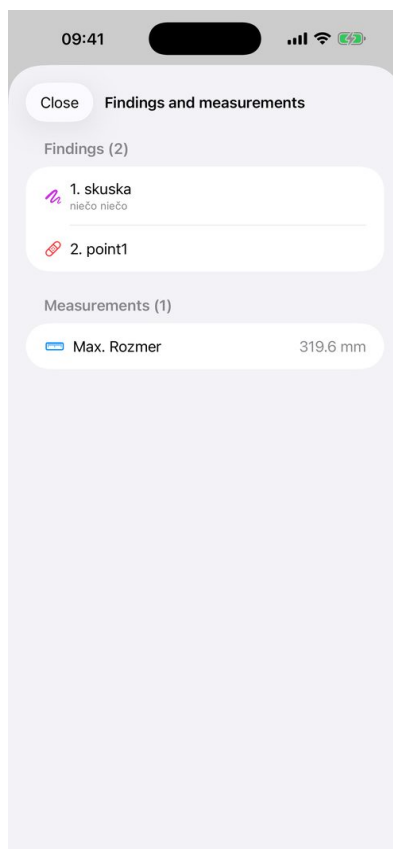

List of findings and saved measurements



Canister — Object / item profile, Detail extent


8. Measuring with the smallest error

Switch the panel to Measure and choose the kind: Distance (2 points, straight line) · Path (along the surface — wound length on a curved body, circumference) · Area (polygon ≥ 3 points — bruise, abrasion) · Volume — box (bounding box from points) · Volume — model (whole mesh, indicative).

- 1 Reference first, everything else after. Measure the 100 mm scale of the target with Distance (line to line, zoomed in), tap Reference, enter 100. The deviation is stored in the report; it should be within ± 1 %. If larger, repeat the scan with the target closer to the subject.
- 2 Zoom in so both ends of the measurement span the screen. Tap edges and clearly defined points (end of the wound, margin of the bruise), not “somewhere near them”.
- 3 Use Distance for dimensions across a flat or slightly convex surface (body length: crown–heel). If “along the surface” also appears under the value, the surface bulges between the points — record both values and say which is which.
- 4 Use Path for wound length on a curved surface and for circumferences: tap points along the line every 2–5 cm; between taps the length adds up along the real mesh surface. A path is never shorter than the straight line.
- 5 Tap Area around the outline of the finding with 6–12 points in order (not criss-cross). Area is a planar projection — on a strongly curved surface (shoulder, calf) split it into two smaller ones and add them.
- 6 Volume is indicative: the box from points gives an upper bound, the model volume depends on mesh closure and surroundings. Report it as indicative.

- 7 Save right after every measurement; the label may stay empty (“Distance 3” is filled in) or match the finding (“wound length no. 1”).
- 8 Repeat key measurements twice from another view; if they differ by more than 2 mm on a detail or 1 cm on a whole body, take the decisive measurement on a scan with extent Detail.
- 9 Units (☐ → Settings): mm, cm, dm, m, inches, feet; the conversion is exact, the package always stores mm.


A measurement is defensible when: (1) the reference is verified with a deviation within 1 %, (2) the line “Scale vs. LiDAR” in the scan is ≈ 1.000 , (3) points were tapped zoomed in on clear edges, (4) on a curved surface you used Path, not Distance, (5) the record states: “aPrehliadky 2.1 3D model, reference verified, deviation x %”. A 0.1 mm display resolution is not accuracy — the reference determines accuracy.

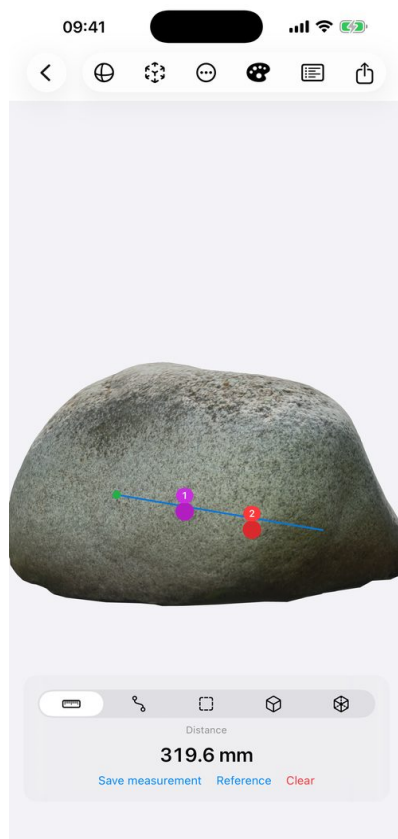

Distance: two points on the surface



Area: four corners of the label (example)

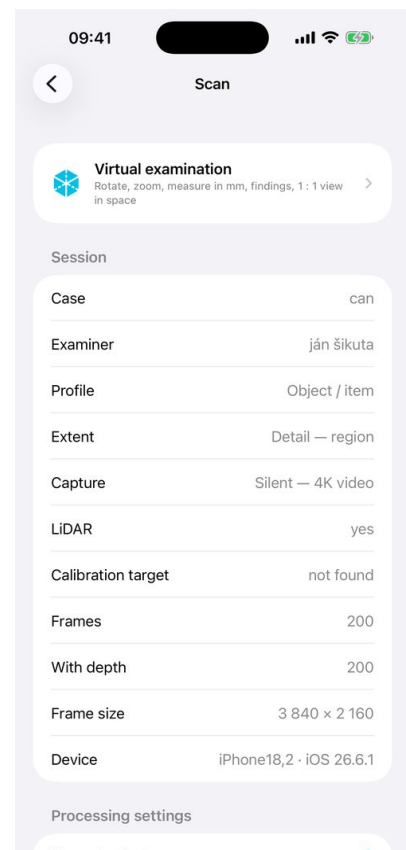


Scan screen: parameters, processing, diagnostics

9. Surroundings from LiDAR — room mesh

With LiDAR the app also collects the ARKit room mesh during capture. After capture a second view appears in the scan, Surroundings from LiDAR: a metric, untextured model coloured by classification (floor, walls, ceiling, table, seat, window, door). It is a direct LiDAR measurement, not photogrammetry and nothing is filled in.

- Suited to the position of the body in the room and distances from walls, doors and furniture; it has its own findings and measurements (different coordinates from the body model) and its own report section with SHA-256.

- Mesh resolution is in the order of centimetres — measure small wounds on the photogrammetric model, the distance of the body from a wall on the mesh.

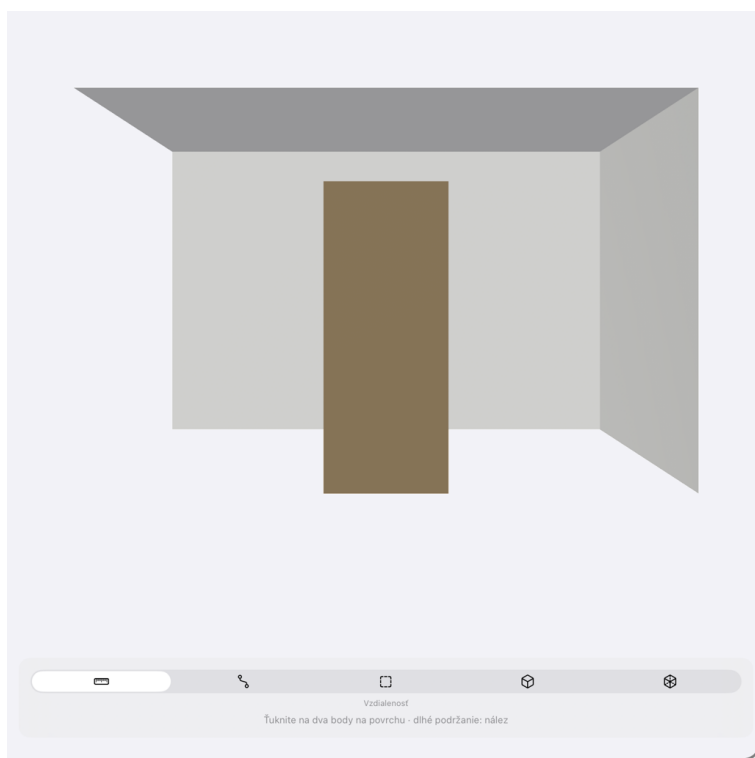

Surroundings from LiDAR: colours by classification (walls, ceiling, door), measurement in metres.

10. Scale — when a measurement is binding

Photogrammetry alone does not know scale. aPrehliadky obtains it from three independent sources and states in the report which one applies:

- 1 LiDAR depth: every frame carries a depth map embedded in the HEIC, the model is built directly in metres. Report line: “HEIC carries embedded LiDAR depth”.
- 2 Camera trajectory (ARKit): camera positions are metric even without LiDAR. After processing they are matched with the photogrammetry poses; the line “Scale vs. LiDAR” should be ≈ 1.000 .
- 3 Reference / calibration target: the only check a court or a colleague can verify by eye on the model. The deviation in % goes into the report.

With the target: scale verified twice — measurements binding. Without the target, with LiDAR: reliable but without independent proof — state “no reference” in the report. Without the target and without LiDAR: indicative; add a reference afterwards by measuring any object of known length in the model.

11. Report, export, archive, Mac

- PDF report: case, examiner, capture parameters, method statement (no NeRF/3DGS/AI), scale verification, findings with positions, measurements in the chosen units, surroundings from LiDAR, integrity, SHA-256 of model and mesh, limitations of the method.
- .f3d package (Forensika 3D; older .aprehliadka files import the same way): manifest, model, findings, certificate, room mesh, processing log, SHA-256 of every file; optionally the frames (allow recomputation elsewhere). The recipient sees on import whether the package was modified. With

the Pro layer the package is additionally signed with the device key (Secure Enclave): the report states “Package signature: VALID, key ..., signed by ...” and the institution name in the header (□ → Pro).

- 3D documentation certificate: Export → Issue certificate (for a finished model). It states the documentation class A–D (A = scale substantiated by two independent routes within 1 %, B = metric without a second proof, C = indicative, D = unverified), the capture trace (camera path, coverage, time), scale verification (LiDAR, trajectory fit, reference, T2 target), the SHA-256 of the model and findings, a signature (Pro) and a QR code. It is also stored in the package (certificate.json) and import verifies it. Methodology: web → Certification methodology (PDF).
- Archive: full package into Files → aPrehliadky → Exporty → Archív; the package SHA-256 is written into the manifest and the report. Frames can then be deleted from the phone.
- Save to... on export opens the system folder picker (iCloud Drive, external drive); on the iPhone exports are also in Files → On My iPhone → aPrehliadky → Exporty.
- Mac (Apple Silicon): double-click a package from the iPhone — examination, measuring, findings and the report work; capture and processing stay on iPhone/iPad.

12. ePrehliadky (ÚDZS) and legal notice

aPrehliadky is designed as a possible add-on to ePrehliadky — the software of the Slovak Health Care Surveillance Authority for electronic post-mortem examinations. The .f3d package, the PDF report and the certificate are attached to the record; the link `aprehliadky://scan/<id>` leads back to the scan. The package format is open (docs/FORMAT.md) and usable in any other system.

aPrehliadky is an independent app by the author (© 2026 Ján Šikuta, MD, PhD). It is not a product of ÚDZS, and ÚDZS neither publishes nor endorses it.

13. Troubleshooting

- The device cannot compute the model: capture, Export → Model, findings and frames, send the package by AirDrop to a supported device, import and start reconstruction.
- I cannot see the model: Scan → Diagnostics → Open in the system viewer. If the model exists it always shows. If not, send the processing log (it is in the package) to jansikuta@me.com.
- Target not found: print at 100 %, target flat, in as many frames as possible, no glare.
- Measurement disagrees with a ruler: verify the reference; if the deviation exceeds 1 %, repeat the scan with the target closer to the subject and a slower walk.
- The phone heats up: normal during processing; leave it on the charger and undisturbed.
- Language: □ → Settings → App language (takes effect after restart).

14. Error codes — meaning and remedy

Every error in the app is shown with a code E-xx, a short meaning and a remedy; the original system message (domain and number) is in square brackets and should be quoted when reporting. The most common one is E30 (“system code 1” from photogrammetry).

- E01 Device does not support capture (ARKit). — Capture on an iPhone/iPad with iOS 18; on a Mac and in the Simulator only import works.

- E02 Frames or depth could not be written. — Free up storage, repeat the capture; if it recurs, restart the device.
- E03 Session without frames. — Capture at least 20 frames before finishing.
- E04 Device cannot compute the model. — Recompute the package with frames on an iPhone/iPad with A17 Pro / M1 or newer.
- E05 The .f3d package is damaged. — Request the package again (AirDrop, not a Mail preview); do not unzip and re-zip.
- E10 Processing cancelled. — Tap Recompute.
- E20 Storage full · E21 File not found (frames deleted after archiving) · E22 No access to the file (open via Files → Share → aPrehliadky).
- E30 Photogrammetry could not process the frames (system code 1 — frames could not be matched/aligned). — 1) Recompute with sensitivity High and extent Whole scene; 2) input Folder without depth (diagnostics); 3) repeat the scan more slowly, with more overlap (three heights, 0.8–1.2 m), no subject movement, no glossy or plain surfaces. Causes: blurred frames, little texture, subject movement, changing light, frames against a window.
- E31 Photogrammetry rejected the input frames. — Try Folder without depth; if it passes, the depth data is damaged (capture again in silent mode).
- E32 Output cannot be written · E33 Not enough space for processing (states how much is needed). — Free up storage: archive old scans and delete their frames.
- E34 Photogrammetry input is not valid (code 2). — Frame diagnostics, recompute with Sample sequence.
- E35 Processing interrupted by the system (code 3/4 — memory/temperature). — Cool down on the charger, close other apps, extent Subject.
- E39 Other photogrammetry failure (system code given). — Sensitivity High; if it recurs, send the processing log.
- E40 Unsupported ARKit configuration (Silent mode, disable surroundings) · E41 Camera/LiDAR unavailable (close other camera apps) · E42 Camera not allowed (Settings → aPrehliadky → Camera) · E43 Tracking failed (stop, textured surface, more light) · E44 Full-resolution frame failed (Silent mode) · E45 Too little detail in view · E46 ARKit write failed (storage) · E49 Other ARKit error (restart the app).
- E60 App Store (donor/Pro) — connection and sign-in; the fee is not required for any feature.
- E99 Unknown error — the domain and number are in brackets; send the processing log (Scan → Processing log → Copy) to jansikuta@me.com.

How to report an error: the code (E30), the text in brackets (e.g. “com.apple.CorePhotogrammetry 1”), iPhone model and iOS version, app version (□), and the processing log. Do not send frames anywhere — the log is enough.

15. Privacy and donor

The app has no network layer. Frames, models, findings and reports stay on the device until you export them yourself. No analytics, no advertising.

aPrehliadky is free and fully functional with no limits — capture, measuring, findings and the report are and will remain free. Anyone who wishes can become a donor (☐ → Donor). The Pro layer (package signing, institution on the report) is optional and is recognised automatically for donors.